

HYGIENIC FILTER HOUSINGS

ISH10i

(Steam and gases)

DESCRIPTION

The ADCAPure ISH10i high efficiency filters are used to remove contaminated particles from gases such as steam and compressed air.

The optimized construction of these units offers low differential pressure at high flow rates.

All sizes are built in two halves and are joined by a sanitary clamp ferrule according to DIN 32676 Series A.

MAIN FEATURES

Clamped body.

Fabricated from AISI 316L stainless steel.

Body machined from bar stock.

Optimized flow paths for low pressure drop and high throughput.

FDA / USP Class VI compliant seals.

STANDARD SURFACE FINISH

Internal wetted parts: $\leq 0,51 \mu\text{m Ra} - \text{SF1}$.

External: $\leq 0,76 \mu\text{m Ra} - \text{SF3}$.

Other surface conditions see TIS.GIA – General information ADCAPure.

Ultrasonic cleaning.

OPTIONS: Air vent and bottom drain connections.
Other connections and dimensions.
Vertical installation design (ISV10i).

USE: Steam, compressed air and other gases.

AVAILABLE MODELS: ISH10i.

FILTER ELEMENTS: Consult respective IS information sheets.

SIZES: 1/2" to 2"; DN 10 to DN 50.




CONNECTIONS: ASME BPE, DIN and ISO clamp ferrules.
Others on request.

INSTALLATION: Horizontal installation with the drain connection pointing downwards.
See IMI – Installation and maintenance instructions.



ISV10i
Vertical inlet/outlet

CE MARKING – GROUP 2 (PED – EUROPEAN DIRECTIVE)	
Ps 16 bar	CATEGORY
1/4" to 2"L – DN 10 to 50L	SEP
2"H – DN 50H	1 (CE marked)

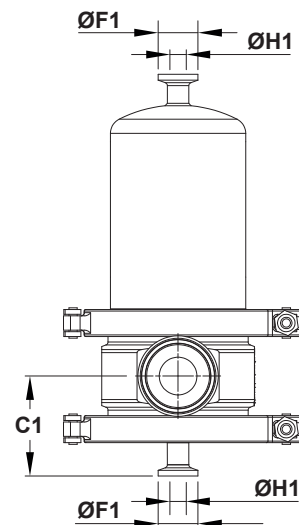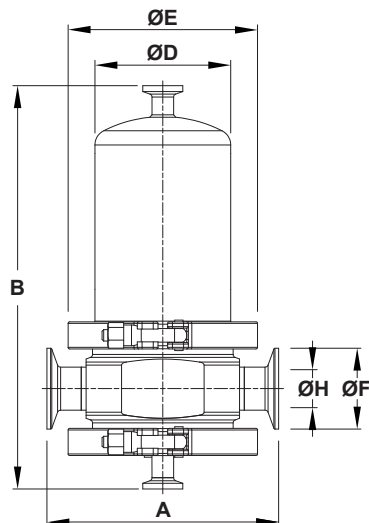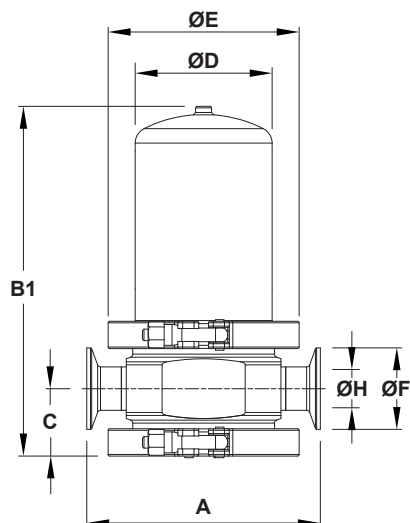

DIMENSIONS – ASME BPE (mm)

SIZE *	A	B	B1	C	C1	ØD	ØE	ØF	ØF1	ØH	ØH1	ISFE SIZE	ISFE QTY.	VOL. (L)	WGT. (kg)
1/2"	140	204	154	34	55	70	104	25	25	9,4	9,4	0310	1	0,31	3,3
3/4"	140	237	187	37	58	70	104	25	25	15,75	9,4	0420	1	0,37	3,6
1"	159	272	222	42	63	85	119	50,4	25	22,1	9,4	0520	1	0,84	5,3
1 1/2"	161	339	288	49	70	85	119	50,4	25	34,8	9,4	0725	1	1,22	6,3
2"L	174	418	369	52	73	104	134	63,9	25	47,5	9,4	1030	1	2,15	7,9
2"H	174	545	496	52	73	104	134	63,9	25	47,5	9,4	1530	1	3,56	8,7

DIMENSIONS – DIN (mm)

SIZE *	A	B	B1	C	C1	ØD	ØE	ØF	ØF1	ØH	ØH1	ISFE SIZE	ISFE QTY.	VOL. (L)	WGT. (kg)
DN 10	140	204	154	34	55	70	104	34	34	10	10	0310	1	0,31	3,4
DN 15	140	237	187	37	58	70	104	34	34	16	10	0420	1	0,37	3,7
DN 20	159	272	222	42	63	85	119	34	34	20	10	0520	1	0,8	5,2
DN 25	174	272	222	42	63	85	119	50,5	34	26	10	0525	1	0,81	5,2
DN 32	176	344	295	49	70	85	119	50,5	34	32	10	0725	1	1,19	6,4
DN 40	189	344	295	52	73	104	134	50,5	34	38	10	0730	1	1,64	7,6
DN 50L	189	418	369	52	73	104	134	64	34	50	10	1030	1	2,32	7,8
DN 50H	189	545	496	52	73	104	134	64	34	50	10	1530	1	3,64	8,6

Remark: Clamp ferrules according to DIN 32676-A (for pipes DIN 11866-A – DIN 11850-2).

DIMENSIONS – ISO (mm)

SIZE *	A	B	B1	C	C1	ØD	ØE	ØF	ØF1	ØH	ØH1	ISFE SIZE	ISFE QTY.	VOL. (L)	WGT. (kg)
DN 08	140	204	155	34	55	70	104	25	25	10,3	10,3	0310	1	0,35	3,3
DN 10	140	237	187	34	55	70	104	25	25	14	10,3	0410	1	0,45	3,5
DN 15	140	242	192	37	58	70	104	50,5	25	18,1	10,3	0420	1	0,46	3,7
DN 20	145	272	222	42	63	85	119	50,5	25	23,7	10,3	0520	1	0,85	5,1
DN 25	145	282	232	42	63	85	119	50,5	25	29,7	10,3	0525	1	0,89	5,1
DN 32	147	344	294	49	70	85	119	64	25	38,4	10,3	0725	1	1,26	6,3
DN 40	160	344	295	52	73	104	134	64	25	44,3	10,3	0730	1	1,95	7,6
DN 50L	173	433	384	57	78	104	134	77,5	25	56,3	10,3	1030	1	2,69	8,4
DN 50H	173	560	511	57	78	104	134	77,5	25	56,3	10,3	1530	1	3,71	9,2

* Suffix L corresponds to low capacity design; suffix H corresponds to high capacity design.

Remarks: Clamp ferrules according to DIN 32676-B (for pipes DIN 11866-B – ISO 1127).



Other sizes on request.


LIMITING CONDITIONS *

Maximum allowable pressure	16 bar
Maximum allowable temperature	200 °C
Minimum allowable temperature	-20 °C
Maximum cold hydraulic test pressure	28 bar
Maximum differential pressure	5 bar

* Other limits on request. Maximum operating conditions may be limited by the filter housing end connections due to normative restrictions.

OPTIONS

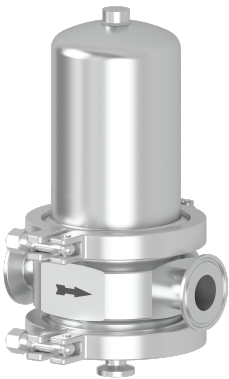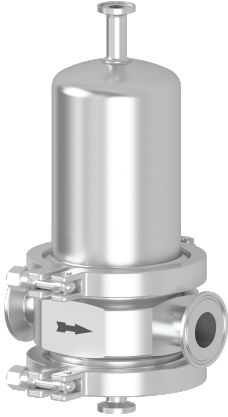BOTTOM DRAIN CONNECTION FOR CLEAN STEAM APPLICATIONS	AIR VENT AND BOTTOM DRAIN CONNECTIONS
	

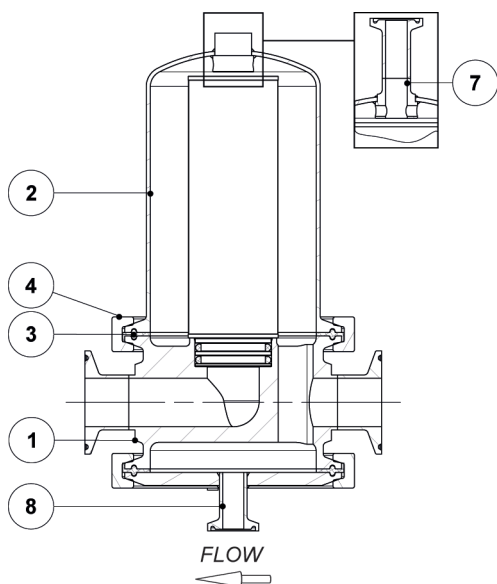

MATERIALS

POS. No.	DESIGNATION	MATERIAL
1	Filter housing	AISI 316L / 1.4404
2	Filter housing cover	AISI 316L / 1.4404
3	* Seal	** PTFE / FPM Envelope
4	Safety clamp	AISI 316 / 1.4401
7	Air vent connection	AISI 316L / 1.4404
8	Drain connection	AISI 316L / 1.4404

* Available spare parts. ** Others on request.

FDA / USP Class VI seals certificate on request.

All filters have a serial number. In case of non-standard filters this number must be supplied if spare parts are ordered.

ORDERING CODES – ISH10i FILTER HOUSING										
FILTER MODEL	ISH10	X	T	X	X	DX	0	008X		
ISH10i – Hygienic steam filter, filter housing	ISH10									
MATERIAL										
AISI 316L / 1.4404		X								
BODY SEALING										
PTFE / FPM			T							
SURFACE FINISH (a)										
Standard surface finish					X					
Mirror mechanical polished external surfaces (SF1)					P					
Electropolished internal wetted parts (SF5)					E					
Electropolished internal wetted parts (SF4)					Q					
SPECIAL FEATURES										
None						X				
Degreased for oxygen						O				
PIPE CONNECTIONS										
Clamp ferrule ASME BPE							DX			
Clamp ferrule DIN (DIN 32676-A)							FX			
Clamp ferrule ISO (DIN 32676-B)							EX			
DRAIN AND VENT CONNECTIONS										
Without vent or drain connection								0		
With bottom drain connection								1		
With bottom drain and vent connection								2		
SIZE										
1/4" or DN 08									008X	
3/8" or DN 10									010X	
1/2" or DN 15									015X	
...									...	
2"L or DN 50L									050L	
2"H or DN 50H									050H	
SPECIAL CONSTRUCTION / ADDITIONAL OPTIONS										
Full description or additional codes have to be added in case of non-standard combination										

(a) Consult TIS.GIA – General information ADCAPure – for further details and other surface finish options.

HOW TO ORDER

1 ADCAPure filter housing ISH10i with PTFE seals and clamp ferrules ASME BPE 2"L – Code: ISH10XTXXDX0050L
 1 ISFE filter element with 5 micron retention rate and EPM seals for the above mentioned filter – Code: ISFE.P1030.05
 Remark: We recommend a second filter element set as a spare part to ensure minimum downtime when replacing the one in use after saturation.